

VIA Anti-Gravity Treadmill™

The AlterG® Via Anti-Gravity Treadmill changes what's possible in rehabilitation and training. With unmatched precision and control, the Via combines revolutionary Differential Air Pressure (DAP) technology, an intuitive touchscreen interface, real-time gait analytics, and live video monitoring – all to enable the ultimate in pain-free and fall-safe body weight supported walking, running, and closed-kinetic chain exercise. Now you can help diverse ortho and neuro patients, active seniors, and athletes to gain mobility, recover smarter, and achieve performance goals sooner.

NEW DESIGN

Accessible Design

Snap-up, magnetized air chamber closure lays flat for ease of entry, whilst a new window design provides critical visual accessibility in a robust, functional design.

Smooth Setup

Open-back cockpit design with light and easy height adjustment leverages inseam set-point measurements and smooth-set locking mechanism to make entry and setup simple.


Easy In, Easy Out

New, low step-up height ensures the easiest possible in-and-out for users.



Enhanced Data

AlterG Stride Smart provides valuable insight into gait and session data, allowing user-progress evaluation with essential data metrics.

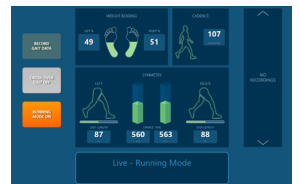

Functional Safety

AlterG's patented DAP technology reduces gravitational forces and provides a fall-safe environment, whilst supporting unimpeded gait mechanics and balance.

VIA X
Anti-Gravity Treadmill

ANTI-GRAVITY PLATFORM



SPECIFICATIONS

| | | | |
|-------------------------|--|---|-------------------|
| Body Weight Adjustment | 100% down to as little as 20% in 1% increments | | |
| Display | 15" | 15" | 21" |
| Treadmill | Flat belt | Flat belt | Slat belt |
| Speed | 0 mph to 12 mph | -3 mph to 12 mph | -10 mph to 18 mph |
| Starting Speed | 0.1 mph | 0.1 mph | 0.1 mph |
| X Package (High Speed)* | Optional (0 mph -15 mph, 0.5 mph Starting Speed, Fit X Graphic) | Optional (-5 mph to 15 mph, 0.5 mph Starting Speed, Via X Graphic) | - |
| Reverse (-3 mph) | - | Included | Included |
| Incline (0% - 15%) | - | Included | Included |
| Live Video Monitoring* | - | Included | Included |
| Stride Smart | - | Included | Included |
| AlterG Assistant | - | Included | Included |
| User Weight | 85 lbs to 400 lbs | 85 lbs to 400 lbs | 85 lbs to 400 lbs |

*Service call required for field upgrade(s).

Learn more at alterg.com or call +1.510.270.5900.